

A woman with long dark hair, wearing a grey tank top and grey underwear, is shown from the waist up. She is looking slightly to the right. The background is a bright, slightly blurred outdoor setting. A semi-transparent blue overlay covers the left side of the image, containing white and black text.

**CM**  
Slim

CONTOURING MASTER™

MAKE  
YOUR  
BODY  
SLIM

 TGA APPROVED

CMSlim™ is the most advanced and intensive magnetic muscle stimulator.

Focused magnetic stimulation passes through all skin and fat in order to effectively stimulate the muscle, providing the most intensive continuous contractions which is ideal for muscle growth, and also inducing apoptosis.


CMSlim™ is an electromagnetic muscle stimulator device designed for aesthetic purpose, having two (2) applicators with higher intensity.

It's cutting-edge technology in non-invasive body contouring, as it not only burns fat, but also builds muscle.

This is key to a toned and healthy physique.

# New Technology in Non-invasive Body Slimming



 **TGA APPROVED**

CMSlim™ is a 93/42/EEC MDD certified, TGA approved medical device. Thanks to its patented water cooling system, magnetic output intensity of 7 Tesla it is the most powerful system in its class today unrivalled by the competitors.




Two applicators are placed on the target muscle area, such as abs, thighs or buttocks. The applicators then generate intense magnetic stimulations that cause involuntary muscle contractions.

# Two Applicators

These contractions trigger the release of free fatty acids, which break down fat deposits and increase muscle tone and strength.

This is very similar to the way of doing a workout.


As working out, the stimulation causes to rebuild and repair muscle tissue, resulting in stronger, firmer muscles.

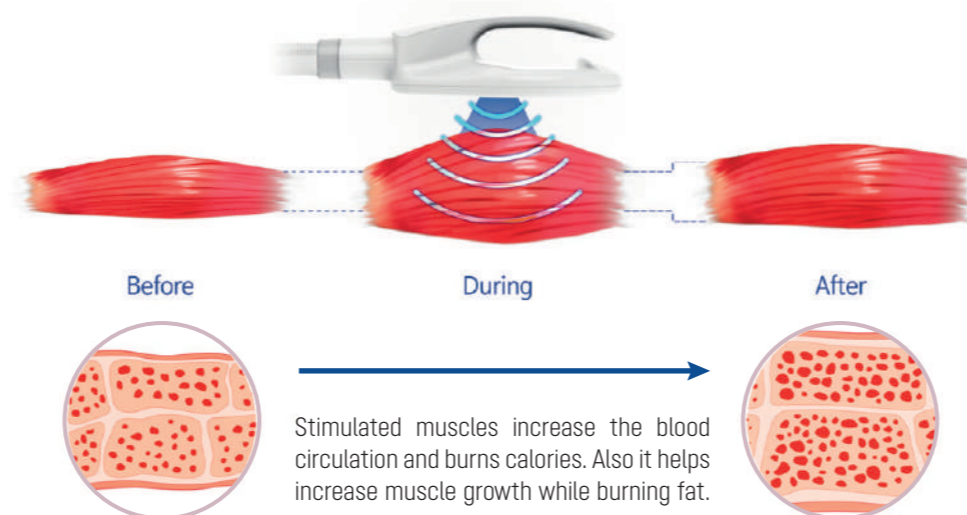




Non-invasive medical technology that is used for strengthening and re-education of muscles via interaction of the magnetic field with the tissue of the patient.

Contrary to voluntary muscle contractions, the supramaximal contractions are independent of brain function. The magnetic stimulation uses a specific range of frequencies that does not allow muscle relaxation between two consecutive stimulations.

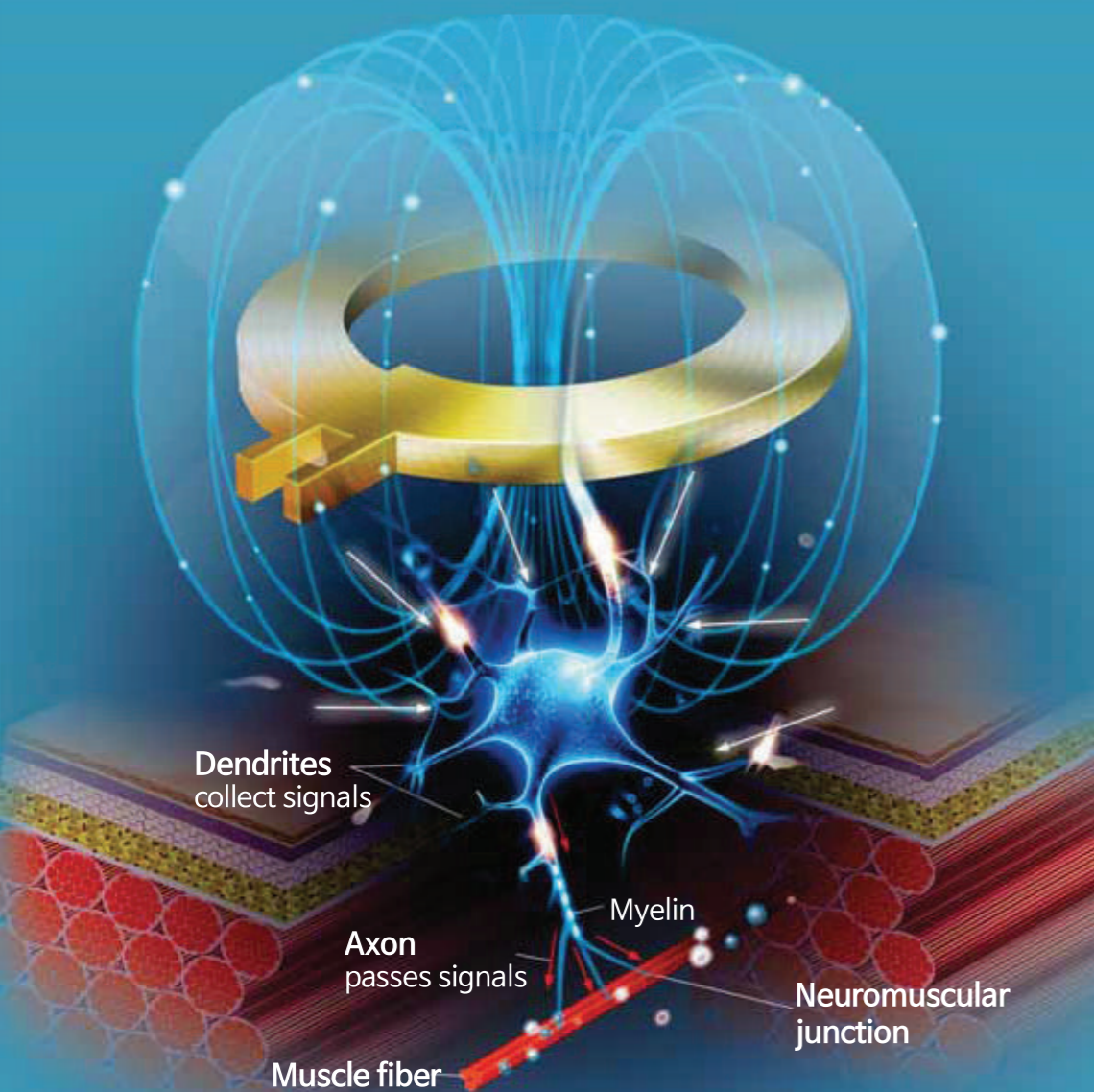

# Proven Technology

The magnetic stimulation uses a specific range of frequencies that does not allow muscle relaxation between two consecutive stimulations. The muscle is forced to remain in a contracted state for multiple seconds. When repeatedly exposed to these high load conditions the muscle tissue is stressed and is forced to adapt.

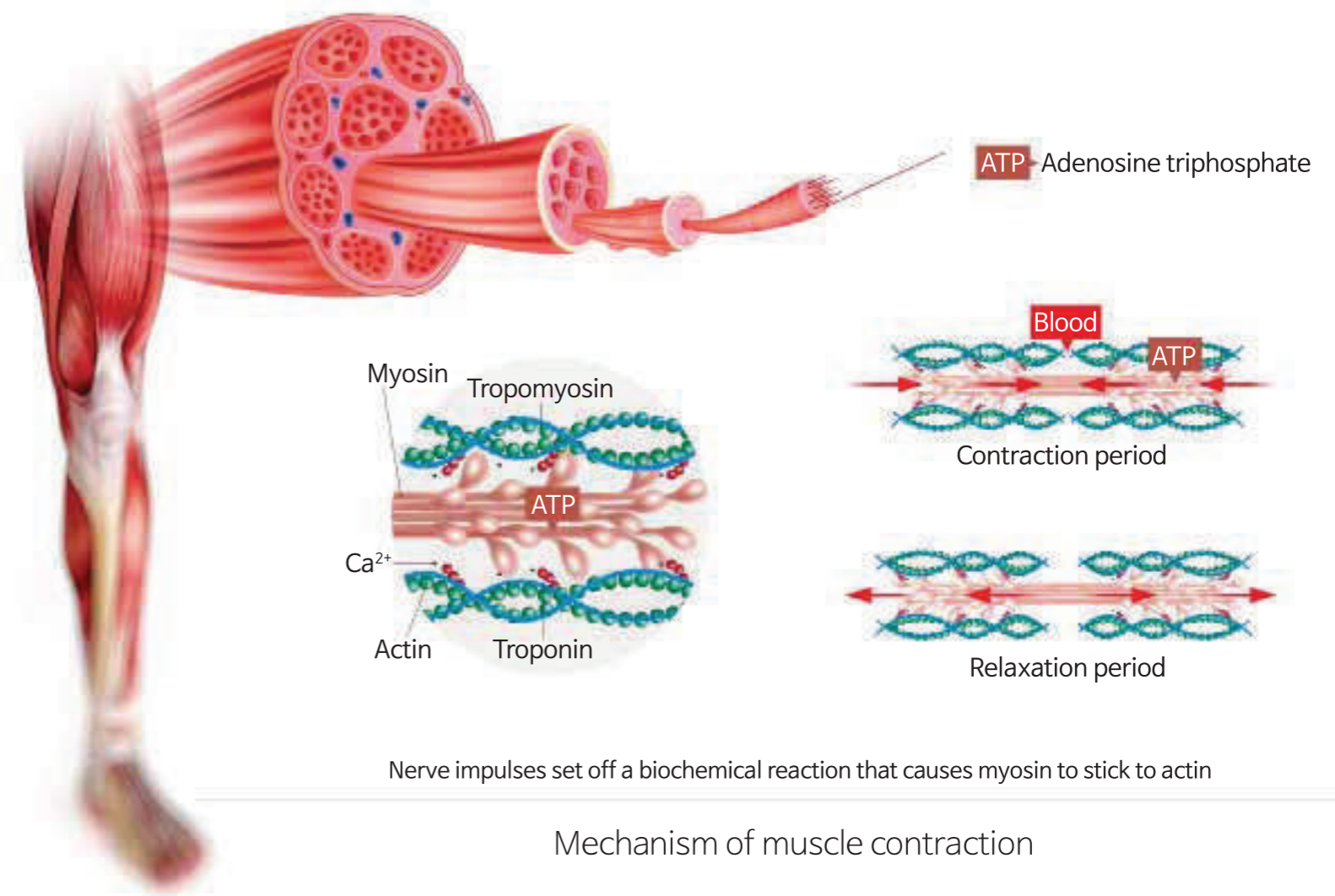

Mechanism of muscle contraction

Recent studies reported that on average 15% - 16% increase in abdominal muscle thickness was observed in treated patients one to two months after magnetic muscle stimulation treatments.






# HI-EMT

## High Intensity Electromagnetic Muscle Trainer

is the medical technology that is used in aesthetic medicine. It uses focused electromagnetic field with a safe intensity level.

Electromagnetic field passes non-invasively through the body and interacts with motor neurons which subsequently trigger supramaximal muscle contractions.

“HI-EMT is used as a non-invasive body contouring treatment for the abdomen, thighs and buttocks.”

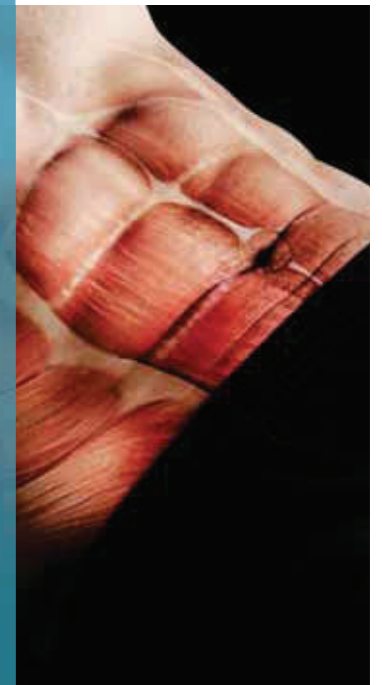

Effect on fat, several recent studies using CT, MRI and ultrasound evaluations have reported approximately 19% reduction in subcutaneous fat layer in patients treated by a magnetic stimulation based device on their abdomen. The principle of cell apoptosis induced by increased concentrations of free fatty acids has been previously observed and demonstrated in numerous research studies.



# Session Information

CMSlim™ treatment session requires no anaesthesia, incisions or discomfort. In fact, patients are able to sit back and relax, while the device performs the equivalent of up to 30,000 painless crunches or squats in 30 minutes.



**Hercules Course**  
Building Muscle

|           |             |        |             |        |             |        |
|-----------|-------------|--------|-------------|--------|-------------|--------|
| Treatment | 1           | 2      | 3           | 4      | 5           | 6      |
| Program   | Hypertrophy | Combo2 | Hypertrophy | Combo2 | Hypertrophy | Combo2 |

**David Course**  
Burning Fat

|           |      |        |          |          |      |        |
|-----------|------|--------|----------|----------|------|--------|
| Treatment | 1    | 2      | 3        | 4        | 5    | 6      |
| Program   | HIIT | Combo1 | Strength | Strength | HIIT | Combo1 |


Single treatment session takes 30 minutes only, and there is no downtime. 2-3 treatments per week, would be enough for perfect result for most people. Generally 4-6 treatments recommended.





# Workout Program

CMSlim™ unique workout program consists of a series of an actual physical training program, which has been ideally designed for the maximum stimulation leading to the best result.

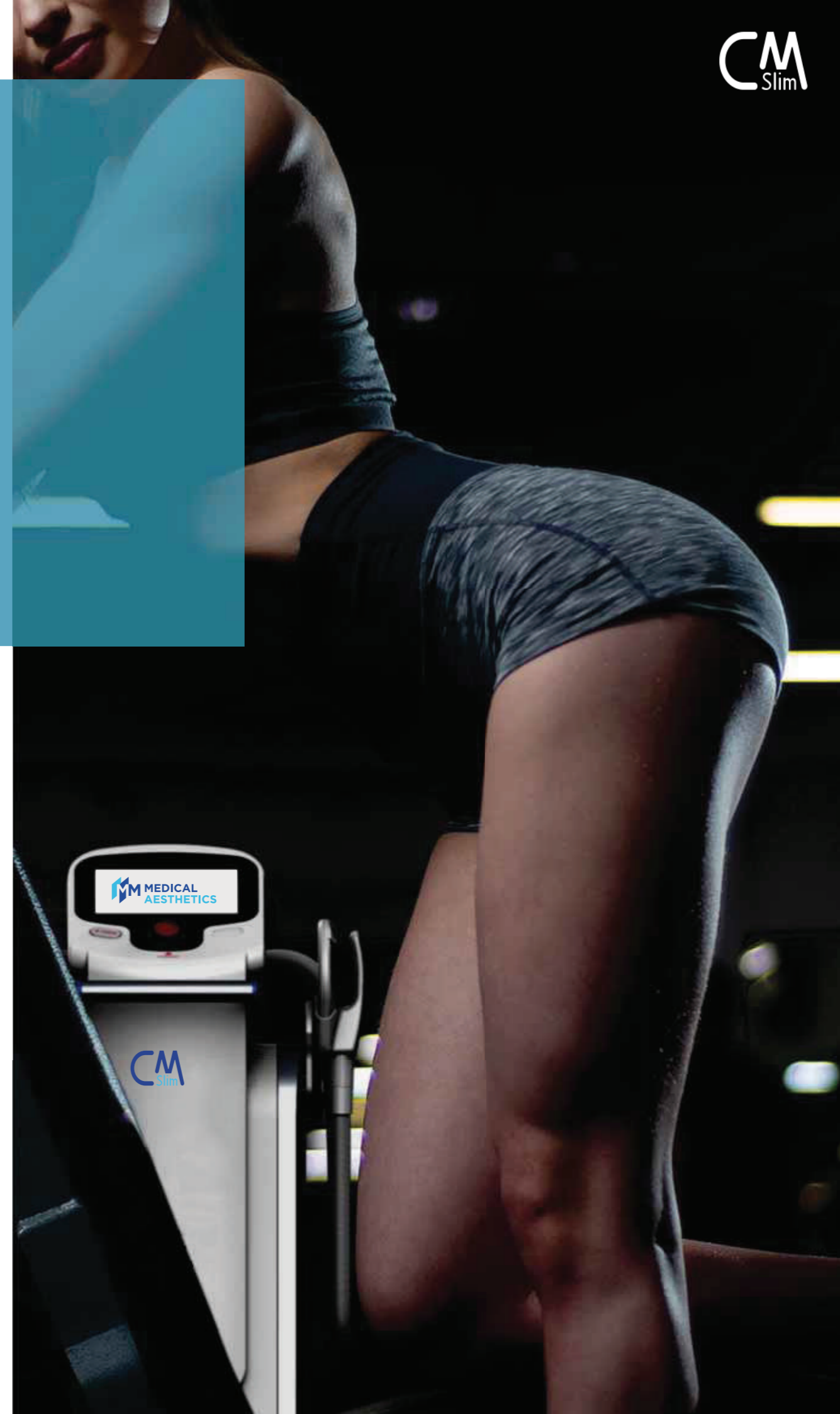

**Venus Course**  
Building Muscle

| Treatment | 1           | 2    | 3      | 4           | 5    | 6      |
|-----------|-------------|------|--------|-------------|------|--------|
| Program   | Hypertrophy | HIIT | Combo2 | Hypertrophy | HIIT | Combo2 |

**Aphrodite Course**  
Burning Fat


| Treatment | 1    | 2       | 3        | 4    | 5      | 6    |
|-----------|------|---------|----------|------|--------|------|
| Program   | HIIT | Combo 1 | Strength | HIIT | Combo2 | HIIT |

Additionally, CMSlim™ is an extremely efficient treatment that takes place over the course of 10-14 days, rather than months, and it also produces more toned abs and uptight buttocks within 2-4 weeks.

If you need more information about workout program, please visit [www.cmslim.com.au](http://www.cmslim.com.au) and send us your enquiry.



# CMSlim™ Advantages

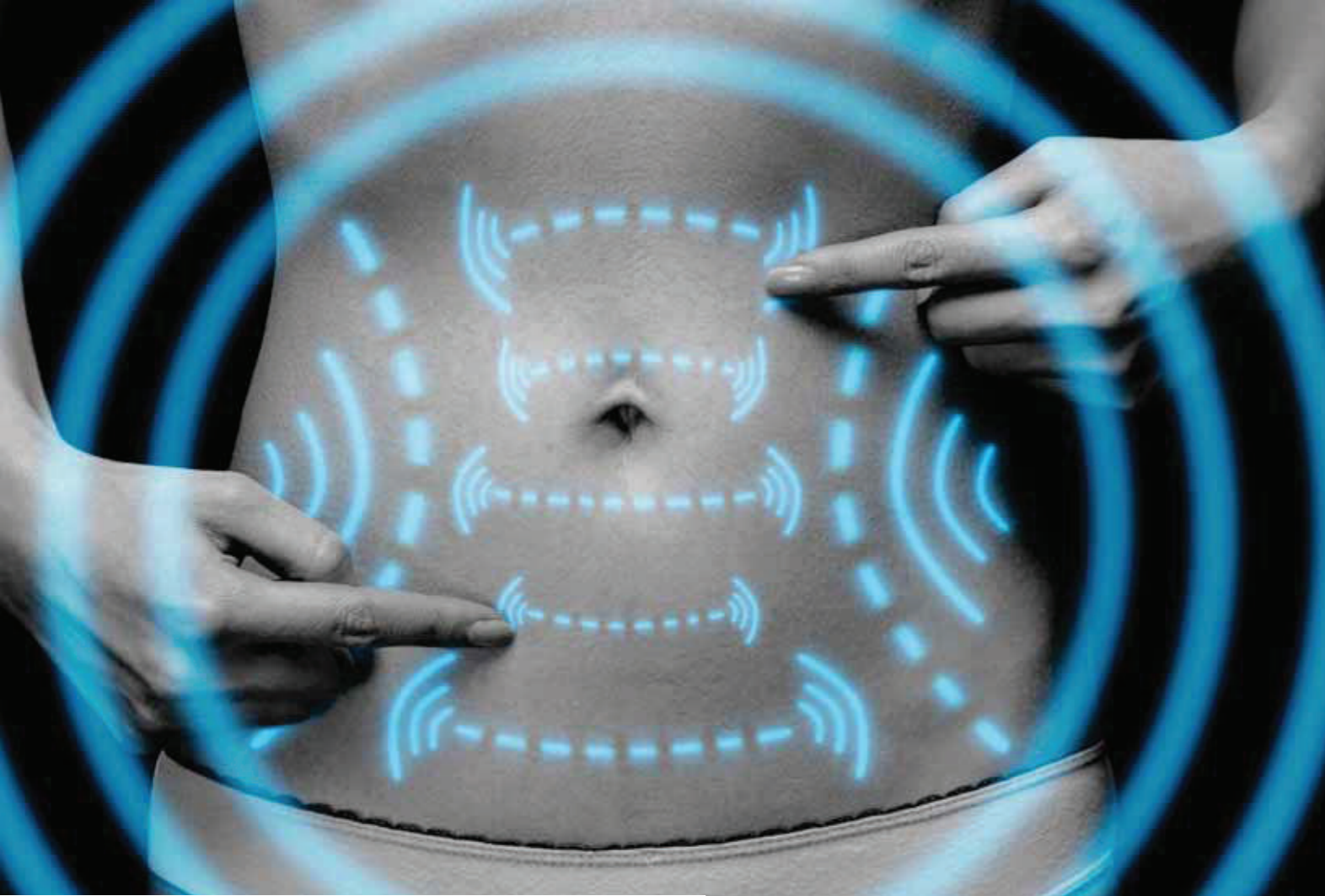

## Output Intensity 7 Tesla


HI-EMT base 7 Tesla output enables that magnetic energy could embrace big skeletal muscles of the human body. That means CMSlim™ muscle re-education process is deep enough like normal muscle contractions. Superficial electrical muscle stimulation is incomparable.



## Cooling Technology

Our optimised water cooling system makes CMSlim™ effective to operate for a very long time without any downtime.

Even further, there is no worry for leaking oils anymore.

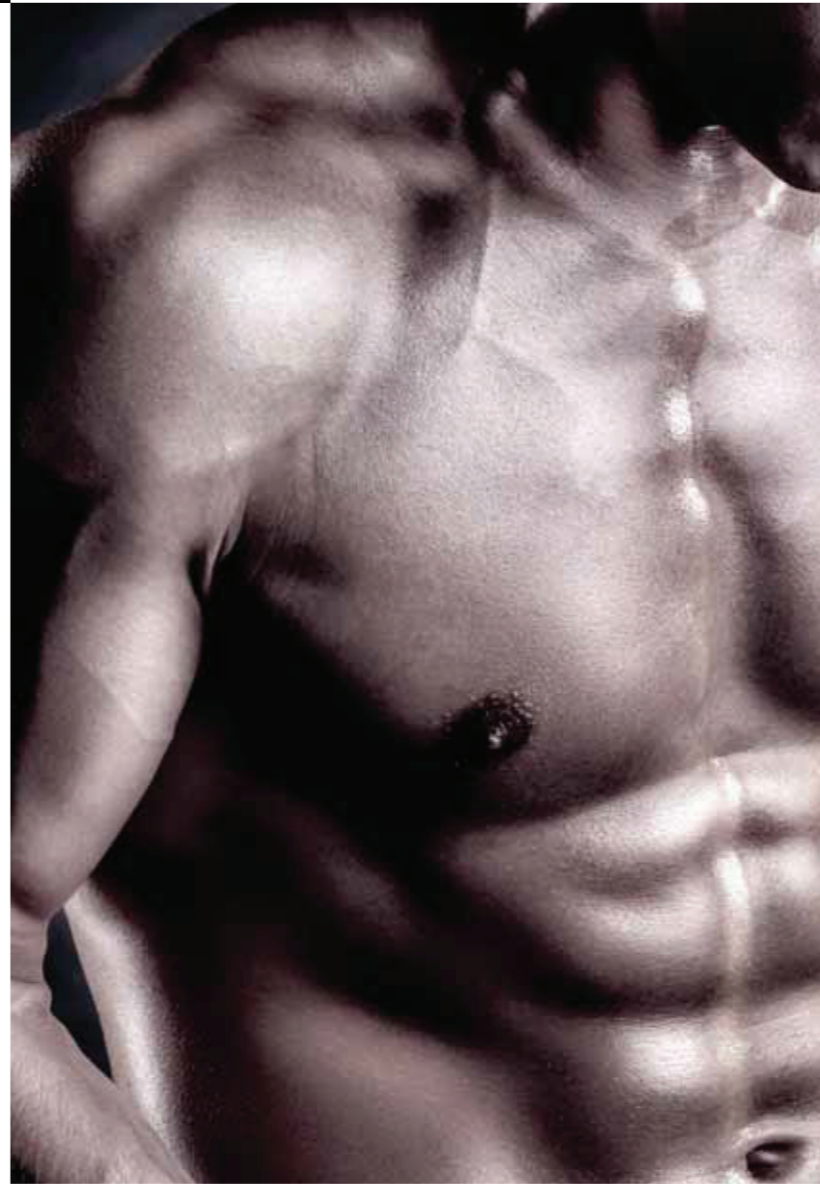

## Manual Training

Our cutting-edge technology makes it possible to adjust magnetic pulse width through to a very wide range. So you don't have to use set presets, you can design special manual treatment programs for your patients' needs.



-  No More Extra Supplies, CMSlim™ has NO CONSUMABLES
-  Perfect Water Cooling System makes CMSlim™ EVEN STRONGER
-  Great Durability – 10,000 Treatment Assurance
-  1 Year Full Coverage Warranty
-  Convenient Operation – EASY & EFFECTIVE PRESET / User Interface

# Maintenance




# Before & After

Even if you aren't in optimal shape, you can experience benefits from CMSlim™ for years to come with regular exercise and balanced eating.



M1\_Before\_Front

M1\_After\_Front

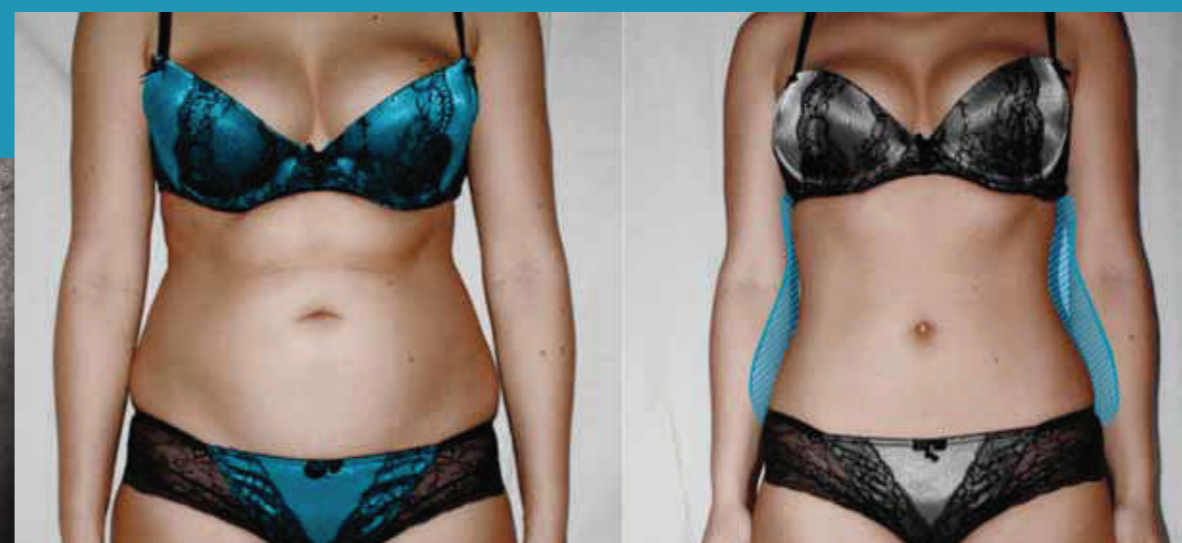

F1\_Before\_Front

F1\_After\_Front



M1\_Before\_Side

M1\_After\_Side



M2\_Before\_Front

M2\_After\_Front



# FAQs

Q A candidate?

A **Anybody**

Q A downtime?

A **Not at all**

Q How long?

A **30 minutes per session**

Q How often?

A **2~3 times a week**

Q How many time?

A **4~6 treatment recommended**

Q How long can I use it?

A **No limit**

Q How long will it last?

A **Extremely long**

Q Additional supplies?

A **No consumables**

Q Warranty?

A **1 year**

Q Manual program?

A **Easy & Effective manual mode**




- ① LCD    ② Emergency stop button    ③ Reset Switch    ④ Applicator (2EA)

# Specification



|                      |  |                       |
|----------------------|--|-----------------------|
| Power                | Power Input                              | 230 V 50/60 Hz        |
|                      | Maximum Power Consumption                | 2.3 kVA               |
| Magnetic Stimulation | Magnetic induction amplitude (Intensity) | 0~7 Tesla             |
|                      | Shape of stimulation pulse               | Biphase               |
|                      | Pulse duration                           | 300 $\mu$ s           |
| Dimension            | Size (height x width x depth)            | 420 x 660 x 1150 (mm) |
|                      | Weight                                   | 70 kg                 |
|                      | Safety type                              | B                     |
| Country of Origin    | Republic of KOREA                        |                       |
| Manufacturer         | DAEYANG MEDICAL CO., LTD                 |                       |



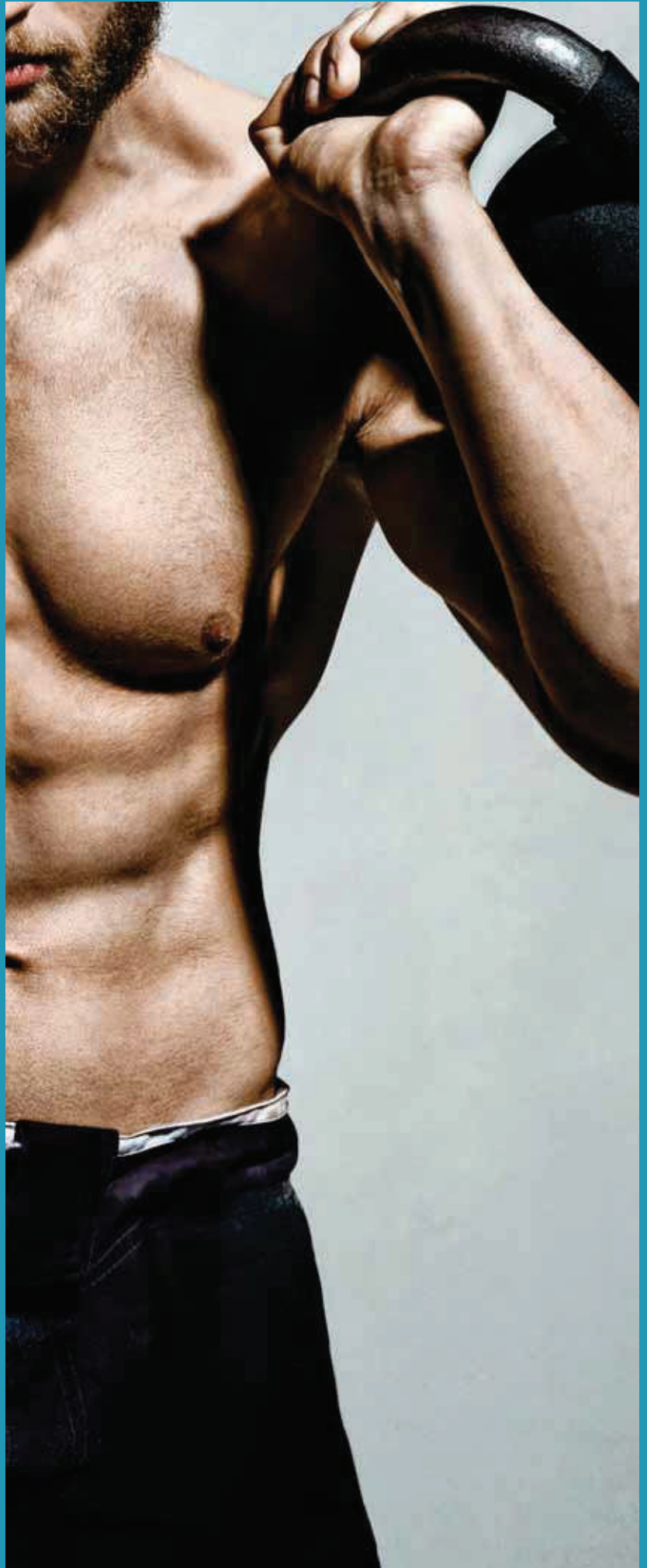

CM  
Slim

Find out more about CMSlim™

 [www.cmslim.com.au](http://www.cmslim.com.au)

 [www.cmslim.com](http://www.cmslim.com)

 [www.mmmedicalaesthetics.com.au](http://www.mmmedicalaesthetics.com.au)

 **MEDICAL  
AESTHETICS**